



August 5, 2026

URGENT FIELD SAFETY NOTICE**iQ200 Series Analyzers and DxU 840m/850m Iris Analyzers**

PRODUCT	REF	Serial Numbers
DxU 840m Iris	C76947	All Serial Numbers
DxU 850m Iris	C49513	
IQ200 SPRINT	700-3322, 700-3322RFB, 700-3325, 700-3326, B48997	
IQ200 SELECT	700-3340, 700-3345, 700-3347, B48998, C10684	
IQ200 ELITE	700-3370, 700-3375, 700-3375RFB, 700-3376, B48999, C10683	

Attention Beckman Coulter Customer,

Beckman Coulter is initiating a field safety corrective action for the product listed above. This letter contains important information that needs your immediate attention. There is no impact on patient results.

ISSUE:	Beckman Coulter has received customer complaints about lamina fluid leaks in iQ200 Series Analyzers and DxU 840m/850m Iris Analyzers. Investigation determined that the leaks are caused by a check valve in the rinse line becoming structurally compromised. This may allow lamina fluid to migrate to internal instrument surfaces and eventually leak onto the table surface and/or floor.
IMPACT:	- Examples of customer observations that may be indicative that the issue is occurring include, but are not limited to: - Visible lamina leakage from the rinse line - Fluid collecting on or around the instrument - iQ Lamina is a colorless, odorless liquid that is not classified as hazardous according to EC 1272/2008 (CLP/GHS), US-OSHA HCS 2012 and UN GHS. Refer to iQ Lamina (PN# 300-7628) Safety Data Sheet for complete safety, handling and disposal information. - Potential health consequences:

Beckman Coulter, Inc.
11800 S.W. 147th Avenue
Miami, FL 33196-9015
USA

www.beckmancoulter.com

FA-26031



	• Slip and fall hazard to laboratory staff from lamina fluid leaking onto the floor, with potential outcomes ranging from bruises or sprains to fractures, head injury, or musculoskeletal trauma. • Discharge of lamina fluid waste in the form of spills, sprays, mist or aerosolization may lead to chemical exposure of operators, service personnel and other users. In a worst case, exposure to <0.15% phenoxyethanol may cause minor eye irritation symptoms, and if they occur, would typically be temporary burning, stinging, tearing, or redness.
ACTION:	• Perform the daily Cleaning Instrument Surfaces maintenance described in Chapter 12, Maintenance Procedure of the DxU 840/850m Instructions for Use (PN# C49320AB), (PN# D38676AA), or iQ200 Instructions for Use (PN# 300-4320CE), (PN# 300-4321EE). During the cleaning, inspect the instrument for signs of lamina leakage, including fluid accumulation on or around the instrument, or evidence of fluid on the laboratory bench or floor. If you determine that there is a leak, discontinue use of the analyzer and clean any spilled fluid in accordance with your laboratory's safety procedure. • Contact Beckman Coulter Technical Support Organization immediately if evidence of leakage is observed. • Continue routine operation, if no leakage is observed, while maintaining increased awareness for the symptoms described in the impact section of this letter. • Beckman Coulter recommends posting this letter on or near iQ200 Series and/ or DxU 840m/850m Iris Analyzers.
RESOLUTION:	• Beckman Coulter is actively investigating this issue. • Beckman Coulter will continue to address any reported lamina fluid leakage through its Technical Support Organization until this issue has been fully resolved.

The national competent authority has been informed of this field safety corrective action.

Please share this information with your laboratory staff and retain this notification as part of your laboratory Quality System documentation. If you have forwarded any of the affected product(s) listed above to another laboratory, please provide them a copy of this letter.

Please complete and return the enclosed Response Form within 10 days so we are assured you have received this important communication.

Beckman Coulter, Inc.
11800 S.W. 147th Avenue
Miami, FL 33196-9015
USA

www.beckmancoulter.com

FA-26031




If you have any questions regarding this notice, please contact our Customer Support Center.

- From our website: <http://www.beckmancoulter.com>

Adverse reactions or quality problems experienced with the use of this product may be reported to the FDA's MedWatch Adverse Event Reporting program either online, by regular mail or by fax.

We apologize for the inconvenience that this caused your laboratory.

Sincerely,

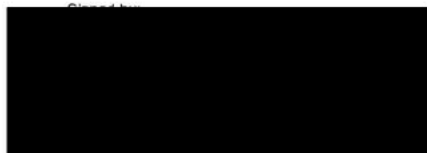

Vice President, U.S. Quality Operations
Beckman Coulter

Enclosure: Response Form

Beckman Coulter, the stylized logo and the Beckman Coulter product and service names mentioned herein are trademarks or registered trademarks of Beckman Coulter, Inc. in the United States and other countries.

Beckman Coulter, Inc.
11800 S.W. 147th Avenue
Miami, FL 33196-9015
USA

www.beckmancoulter.com

FA-26031
